



by RADIANCY™

Facial SkinCare Device With
LHE[®] MicroPhototherapy


Quick Guide

The Quick Guide does not replace the User Manual.
Please read the entire User Manual before treatment.




Know your FSD Plus


Main Unit:

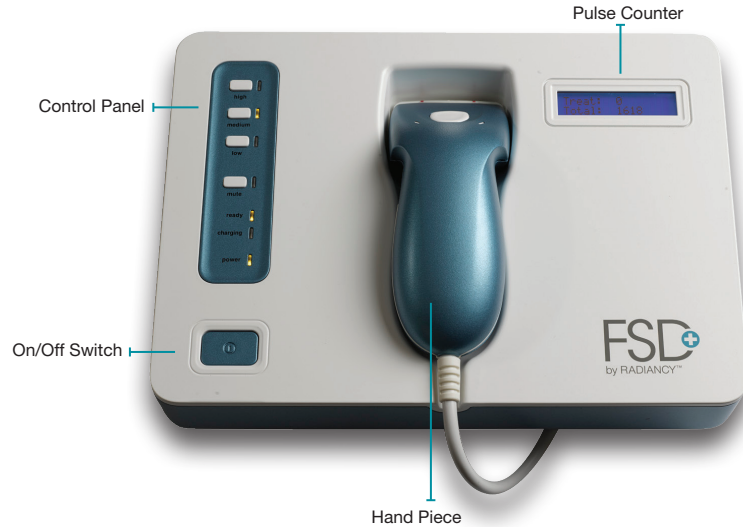


Control Panel:

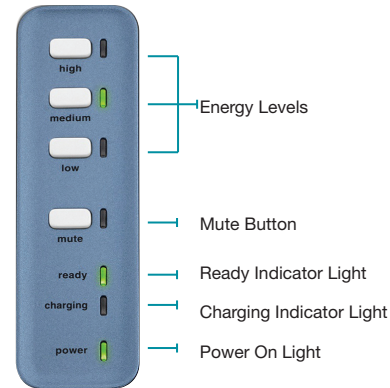

Hand Piece:

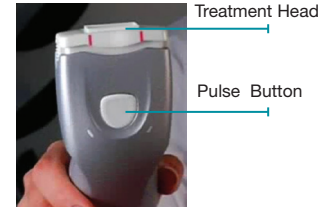

Accessories:

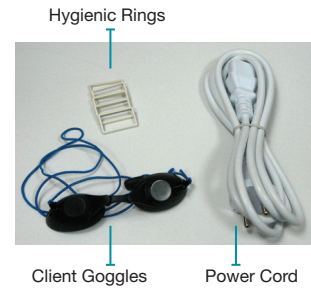

Before You Start

Skin Type Guidelines

- Make sure your client doesn't match any of the exclusion criteria. Refer to the User Manual for further information.
- Determine your client's skin type according to the Fitzpatrick Skin Type Classification below.
- **Do not** treat clients with Skin Type 6.
- Skin Types 4 & 5 should **only** be treated with the low energy setting.

Fitzpatrick Skin Type Classification

					
1	2	3	4	5	6
Always burns, never tans	Always burns, less than average tan	Always burns, average tan	Rarely burns, tans easily	Moderately pigmented, always tans	Deeply pigmented, never tans

Note: Darker skin may be more sensitive to the energy emitted from the lamp and therefore requires a lower energy level setting.

Test Session

Always test your client's reaction to FSD Plus before beginning the actual treatment. Deliver 1-3 pulses in discreet areas and wait at least 24 hours (up to 48 hours for skin types 4-5) for any possible negative reactions.

Use the settings below for the initial test.

Skin Type	Energy Level
1-3	Low-Medium
4-5	Low

Treatment Preparations

1. Make sure the skin is completely clean and dry before beginning treatment.

2. Cover your client eyes with the Goggles.



3. Attach a hygienic ring to the handpiece opening. Use a separate hygienic ring for each client.



4. Attach the power cord to the back of the device and plug the cord into an electrical outlet.



Treatment

1. Turn FSD Plus on.



2. Set the energy level according to the test session. Wait for the "Ready" Indicator light and beep sound. You are now ready to begin treatment.



3. Press the handpiece firmly but gently on the skin. Always make sure:
 - The lamp opening is in full contact with skin.
 - To stretch the skin where you apply the Handpiece.



4. Press and hold the pulse switch down until the handpiece stops flashing or you hear a beep.



5. After the pulse sequence is finished, completely remove the handpiece away from the skin.



6. Apply the next sequence only when the "Ready" Indicator light is back on.



7. Apply the pulses from right to left, in a sequential, horizontal pattern.



8. After covering the entire area, wait one to two minutes to allow the skin to cool.

9. Repeat the procedure: Return to the starting point. This time, hold the handpiece vertically, in the opposite direction of the first pass.



10. When treatment is complete, turn the device off by pressing the power switch.



IMPORTANT NOTE

If your client feels the pulses are too hot, remove the handpiece immediately and wait until the "Ready" indicator light turns on again. If possible, lower the power level. If not, remove after a few pulses.

Notes:

- **Do not** place the handpiece on the skin before the "Ready" Indicator light is on and/or you hear the audible signal.
- Apply pulses with **no more** than a 20% overlap.
- **Stop** if there is full redness of the skin, blistering or burns.
- **Wait** 5 minutes between two consecutive uses.
- **Please read** the "Safety Instructions" section in the User Manual.

After Treatment

Complete the facial for the ultimate affect.

- Apply a cooling, hydrating mask.
- Massage with a non-oily lotion.
- Apply moisturizers (non-stringent, non-alcoholic).
- Apply a moisturizing sun screen (SPF 30 or higher).



RADIANCYTM
Smart Skin Solutions